
AT51400/1/2/3/4/5/6/7/8/9/10/11/12/15

**INSTRUKCJA UŻYTKOWANIA
PIŁKA REHABILITACYJNA
USER MANUAL
REHABILITATION BALL
BENUTZERHANDBUCH
REHABILITATIONSBALL
UŽIVATELSKÁ PŘÍRUČKA
REHABILITAČNÍ MÍČ
NÁVOD NA POUŽITIE
REHABILITAČNÁ LOPTA
MANUEL D'UTILISATION
BALLE DE RÉÉDUCATION
GEBRUIKSAANWIJZING
REVALIDATIEBAL
MANUAL DE USUARIO
PELOTA DE REHABILITACIÓN
MANUALE D'USO
PALLA RIABILITATIVA
BRUKSANVISNING
REHABILITERINGSBOLL**

MD

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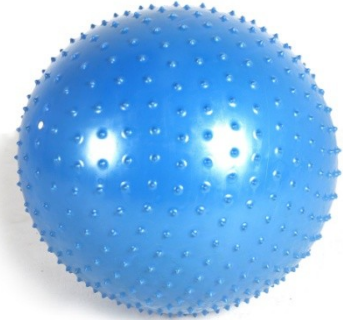EN

Thank you for purchasing our product. Please read the user manual carefully. Do not use the product before familiarizing yourself with this manual and understanding it. If you do not understand the warnings, notes, and recommendations, please contact a healthcare professional or seller to avoid product damage or personal injury.

WARNING: Check all parts for damage that may have occurred during transport. If such damage is noticed, DO NOT use the product. More information is available from the seller.

TARGET PATIENT GROUP

Individuals dealing with diseases, dysfunctions, or injuries for whose treatment/rehabilitation or compensation this product is intended (see the intended use section of this manual). The product may be purchased independently by the user or based on recommendations from a doctor, therapist, or other specialist. Whether purchased independently or recommended by a doctor/therapist/other specialist, available sizes/required features and product variants, indications, and contraindications for use, as well as information provided by the manufacturer, should be considered.

APPLICATION/INDICATIONS

Balls are intended for general fitness and rehabilitation exercises. They are used in fitness exercises, rehabilitation, posture correction, and in prenatal classes. Exercise with a gym ball strengthens muscles responsible for maintaining proper posture.

Using a rehabilitation ball in core-strengthening exercises yields excellent results, helping to maintain correct posture and relieve the spine.

Exercises with a gym ball can be performed anywhere: in an apartment, home gym, outdoors, or in fitness clubs.

Usage:

General fitness exercises

Rehabilitation

Posture correction

Prenatal classes

CONTRAINDICATIONS

Physical or mental limitations (e.g., vision impairment) that prevent safe handling of the product.

Acute injuries or traumas:

People with acute spine, joint, muscle, or other body part injuries should avoid exercises on the ball, as they may worsen the injuries and cause pain.

Spinal diseases:

Patients with serious spinal conditions (e.g., intervertebral disc herniation, scoliosis, osteoporosis) should avoid using the ball without consulting a doctor or physical therapist, as improper exercises may worsen health conditions.

Balance and coordination disorders:

People with balance or coordination issues (e.g., elderly individuals, patients with neurological conditions) should avoid exercises on the ball or perform them only under specialist supervision to minimize the risk of falling.

Acute inflammatory conditions:

In cases of acute joint, muscle, or tendon inflammation, using the ball may worsen symptoms and deteriorate health conditions.

Hypertension and heart disease:

Individuals with uncontrolled hypertension or severe heart disease should avoid exercises on the ball without prior consultation with a doctor, as some positions or exercises may strain the cardiovascular system.

Advanced pregnancy:

Women in advanced pregnancy should avoid certain exercises on the ball that may cause excessive abdominal muscle tension or increase the risk of falls. It is always recommended to consult a gynecologist or midwife before starting exercises on the ball during pregnancy. Dizziness or other neurological disorders:

People who suffer from frequent dizziness, epilepsy, or other neurological diseases should avoid exercises on the ball without specialist supervision.

LIST OF COMPONENTS

Rehabilitation ball, user manual

FEATURES AND DIMENSIONS:

AT51400 "Spiky" massage ball, for rehabilitation - 5CM, PVC

AT51401 "Spiky" massage ball, for rehabilitation - 6CM, PVC

AT51402 "Spiky" massage ball, for rehabilitation - 7CM, PVC

AT51403 "Spiky" massage ball, for rehabilitation - 8CM, PVC

AT51404 "Spiky" massage ball, for rehabilitation - 9CM, PVC

AT51405 "Spiky" massage ball, for rehabilitation - 10CM, PVC

AT51406 Rehabilitation ball with ABS system - 55CM, PVC, MAXIMUM USER WEIGHT 100KG, NON-SLIP

AT51407 Rehabilitation ball with ABS system - 65CM, PVC, MAXIMUM USER WEIGHT 100KG

AT51408 Rehabilitation ball with ABS system - 75CM, PVC, MAXIMUM USER WEIGHT 100KG

AT51409 Rehabilitation ball with ABS system - 85CM, PVC, MAXIMUM USER WEIGHT 100KG

AT51410 Rehabilitation ball - 55CM, PVC, MAXIMUM USER WEIGHT 100KG

AT51411 Rehabilitation ball - 65CM, PVC, MAXIMUM USER WEIGHT 100KG

AT51412 Rehabilitation ball - 75CM, PVC, MAXIMUM USER WEIGHT 100KG

AT51415 "Peanut" rehabilitation ball - 100*50CM, PVC, MAXIMUM USER WEIGHT 100KG

This symbol indicates the maximum user weight.

WARNING:

In the event of a "serious incident" related to the product that directly or indirectly led, could lead, or may lead to any of the following events:

- a) death of the patient, user, or other person or
- b) temporary or permanent deterioration of the health of the patient, user, or other person or
- c) serious public health threat

such a "serious incident" must be reported to the manufacturer and the competent authority of the member state in which the user or patient resides. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices, and Biocidal Products.

WARNING:

In the event of pain, allergic reactions, or other troubling, unclear symptoms associated with the use of the medical device, please consult a healthcare professional.

WARNING: Flammable product.

WARNING: Not intended for use by children under 3 years old.

WARNING: Do not use the product in any way other than its intended use.

CLEANING

The product should be cleaned with a soft cloth and water. Do not use chemicals or strong detergents. Do not wash. Do not iron. Do not bleach.

STORAGE

Store the product in a dry place, away from direct sunlight.

DISPOSAL

Dispose of as household waste in accordance with applicable disposal regulations.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. AN TAR Sp. J., 03-068 Warsaw, ul. Zawiślńska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period AN TAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. AN TAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to AN TAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawiślńska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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